

Kite EMCCD

Digital Monochrome Scientific Interline EMCCD



Actual Size

Key Features and Benefits

Ultra low noise readout for improved image quality!

- **658 x 496 EMCCD sensor.**
Enables optimum image resolution in low light imaging applications
- **B/W EMCCD technology.**
Enables high sensitivity imaging with 1000x on-chip gain
- **16 bit CameraLink output.**
Provides wide dynamic range
- **53% QE from Virtual Phase sensor.**
Optimum Photon collection
- **Frame Interline Transfer (FIT).**
No mechanical shutter required, less smear at shorter exposures
- **High frame rate imaging.**
Enables capture of fast events

Resolution 658 x 496

Readout Noise <1e Noise

Frame Rate 50Hz

CameraLink 16 bit




..... Capturing Tomorrow

Specification for KITE EM247 CL

Sensor	Texas Instruments TC247SPD
Sensor Type	1/2" Frame Interline Transfer (FIT) Impactron
Active Pixel	658 x 496
Pixel Size	10µm x 10µm
Active Area	6.58mm x 4.96mm
Full Well Capacity	24000 electrons
Shift Register Well Depth	100000 electrons
Non Linearity	< 1%
Readout noise	< 1 electrons with EM gain ON < 20 electrons with EM gain OFF
Dynamic Range	65 dB
Frame Rate	50Hz
Dark Current	< 1e / pix / sec
Digital Output Format	16 bit CameraLink (base configuration)
Peak Quantum Efficiency	53% @ 530nm
Spectral response	350 - 1100nm
Cooling	-20°C with ambient air @ +20°C
Binning	1x1, 2x2, 3x3, 4x4, 5x5
Antiblooming protection	Yes
Lens Mount	C mount
Synchronisation	Trigger IN and OUT - TTL compatible
Power Supply	12V DC ±10%
Total power consumption	< 12W
Operating case temperature	-20°C to +55°C
Storage Temperature	-30°C to +85°C
Dimensions	97mm x 68mm x 61mm
Weight (no lens)	< 550g

Quantum Efficiency



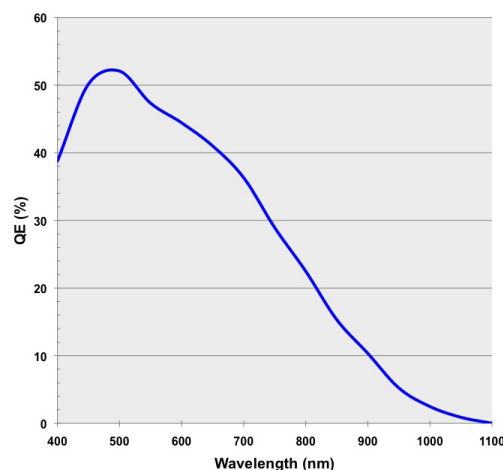


Sample Applications

- Adaptive Optics and Astronomy
- Calcium signaling
- Fluorescence imaging / spectroscopy
- Flow cytometry
- FRET / FRAP / TIRF
- Genome sequencing
- High content screening
- High resolution fluorescence imaging
- Hyperspectral imaging
- LIBS
- Live Cell Imaging
- Single molecule detection
- Solar Cell Inspection
- X-ray tomography

Raptor Photonics Limited reserves the right to change this document at any time without notice and disclaims liability for editorial, pictorial or typographical errors.

Dimensions

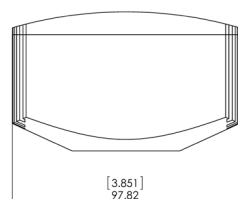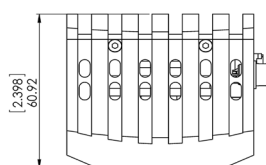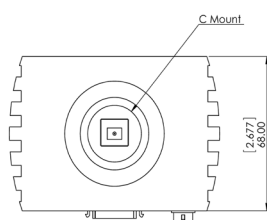

Unit: [inch] mm

Power connector on camera: Hirose HR10-7R-4P

Cable connector (socket): Hirose HR10-7P-4S

Ordering Information

Camera

KITE EM247 digital B/W camera K1247-CL

KITE Power Supply Cable RPL-FA-CBL

Optional Accessories

EPIX(R) base CL card RPL-EPIX-EB1

EPIX(R) base Notebook CL card RPL-EPIX-ECB1-34

EPIX(R) XCAP STD software RPL-XCAP-STD

CameraLink Cable, 2m¹ RPL-CL-CBL-2M

Optical Visible lenses² RPL-xx-xxxx

Note 1: Longer CL cable available

Note 2: Please consult us to check our range of lenses

Equipment may require UK Government authorisation for export purposes

Document #: K1247-CL 0113R1